



BRASELER
USA® SURGICAL

AXIOM™
ULTRASONIC BONE CEMENT REMOVAL SYSTEM

Rapidly Melts, Fragments & Selectively Removes Cement in Hip & Knee Revisions

The Axiom™ Ultrasonic Removal System uses adaptive ultrasonic energy to melt polymethyl methacrylate bone cement (PMMA) at the blade tip, minimizing leverage on bone. The high frequency softens cement allowing selective and rapid cement removal, preserving the surrounding patient anatomy.

KEY BENEFITS

- ✓ Cement removal with minimal force
- ✓ Enhanced procedural guidance
- ✓ Real-time audible & tactile feedback

FEATURES

- ✓ Individual blade optimization stored in console
- ✓ Continuous or pulsed ultrasonic power
- ✓ Tailored power settings from 0-100%





PIERCER

Creates straight relief channels with side vents to evacuate melted cement




SCRAPER

Scoops melted debris, polishes canal walls



ACETABULAR

Curved profile frees cement 360° around cup

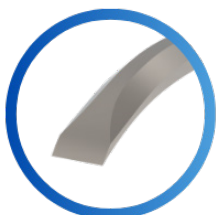

CHISEL

Ultra-thin blade to separate cement mantle



BLADES

Product Number	Item Description	O.D. (mm)	Length (mm)
BRAXM-FP6LH	Full Piercer 6x270, Helical	6	270
BRAXM-FP10S	Full Piercer 10x100	10	100
BRAXM-FP10M	Full Piercer 10x190	10	190
BRAXM-FP10L	Full Piercer 10x300	10	300
BRAXM-HP10S	Half Piercer 10x100	10	100
BRAXM-HP10M	Half Piercer 10x190	10	190
BRAXM-FS6S	Full Scraper 6x100	6	100
BRAXM-FS6M	Full Scraper 6x190	6	190
BRAXM-FS8S	Full Scraper 8x100	8	100
BRAXM-FS8M	Full Scraper 8x190	8	190
BRAXM-HS6S	Half Scraper 6x100	6	100
BRAXM-HS6M	Half Scraper 6x190	6	190
BRAXM-HS8S	Half Scraper 8x100	8	100
BRAXM-HS8M	Half Scraper 8x190	8	190
BRAXM-HS10S	Half Scraper 10x100	10	100
BRAXM-AG6S	Acetabular Gouge 6x100	6	100
BRAXM-AG6L	Acetabular Gouge 6x120	6	120
BRAXM-AG16M	Acetabular Gouge 16x110	16	110
BRAXM-C5R	Chisel 5.5x100 Round	5.5	100
BRAXM-C6R	Chisel 6x100 Round	6	100
BRAXM-C6F	Chisel 6x100 Flat	6	100
BRAXM-EXT	Extender	-	125

SYSTEM COMPONENTS

Product Number	Item Description
BRAXM-UGK	Axiom Ultrasonic Generator Kit
BRAXM-HM	Ultrasonic Handle
BRAXM-UG	Ultrasonic Generator
BRAXM-CWW	Counter Wrench + 8MM
BRAXM-STC	Sterilization Case
BRAXM-FS	Foot Pedal

Call 800.535.6638 to order today!
Visit us online at BrasselerSurgical.com